Mode d'emploi

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
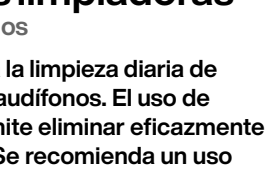
EN Cleaning wipes for hearing aids

Wipes for daily cleaning of all types of hearing aids. Using the wipes effectively removes ear wax. Daily use is recommended.

Precautions for use: External use. Keep out of reach of children. After use, wipe the hearing aids using a clean and absorbent cloth. Wash your hands thoroughly after handling.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.

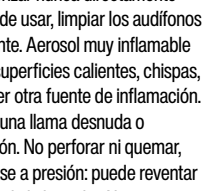EN Gas-free cleaning spray for hearing aids

The spray is intended to remove ear wax from hearing aids. Daily use is recommended.

Precautions for use : External use. Keep out of reach of children. Never spray directly onto the hearing aids. After use, wipe the hearing aids using a clean and absorbent cloth.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.


EN Gas cleaning spray for hearing aids

The spray is intended to remove ear wax from hearing aids. Daily use is recommended.

Precautions for use: Usage externe. Tenir hors de portée des enfants. Ne jamais pulvériser directement sur les aides auditives. Après usage, essuyer les aides auditives avec un tissu propre et absorbant. External use. Keep out of reach of children. Never spray directly onto the hearing aids. After use, wipe the hearing aids using a clean, absorbent cloth. Highly inflammable aerosol - Keep away from heat, hot surfaces, sparks, naked flames and all other sources of ignition. Do not smoke. Do not spray onto a naked flame or any other sources of ignition. Do not pierce or incinerate: even when empty. Pressurised container: may burst under the effects of heat - Protect from direct sunlight. Do not expose to temperatures in excess of 50°C/122°F Causes severe eye irritation. If medical advice is sought, keep the container or label on hand. Wash hands thoroughly after handling. If in contact with the eyes: rinse carefully with water for several minutes.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.

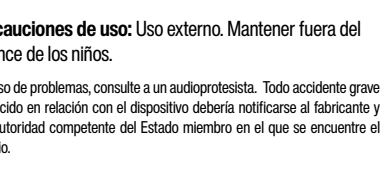
EN Drying Capsules for hearing aids

Capsules to prevent damage to electronic components caused by humidity. Open the battery compartment if the hearing aid has one. Close the cup after each use. Daily use is recommended. Replace the capsule when it is almost white. The service life is about 1 month (drying capsules) / 3 months (drying capsules ultra-powerful).

Precautions for use : Do not pierce the capsule. Do not swallow. Keep out of reach of children.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.

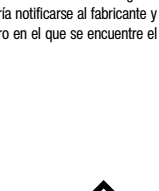EN Cleaning tablets for tubes and earmolds

Cleaning tabs for ear tips. 30 minutes are enough to effectively remove ear wax. Place the cleaning tab in the cup two-thirds full of warm water. After cleaning, rinse the earmold with clean water and dry it. Dispose of the solution after use. Full cleaning of the ear tips is recommended once or twice per week.

Précautions d'emploi : External use. Single use. Keep out of reach of children. Harmful if swallowed. Causes severe skin burns and eye damage. Wash your hands thoroughly after handling. Harmful to aquatic life - Avoid release to the environment. If medical advice is sought, keep the container or label on hand. If swallowed: Call a poison centre / a doctor. Rinse the mouth. Do not induce vomiting. If in contact with the eyes: carefully rinse with water for several minutes.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.

EN Air ball for tubes and vents

The airball is designed to dehumidify tubes and vents.

Precautions for use: External use. Keep out of reach of children.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.

EN Drying spray for hearing aids

The drying spray is indicated for use in dehumidifying hearing aids. Spray using short presses to dry the ear tip. Make sure to remove all residual water from the tube and tip. Daily use is recommended.

Precautions for use: External use. Keep out of reach of children. Highly inflammable aerosol - Keep away from heat, hot surfaces, sparks, naked flames and all other sources of ignition. Do not smoke. Do not spray onto a naked flame or any other sources of ignition. Do not pierce or incinerate: even when empty. Pressurised container: may explode under the effects of heat - Protect from direct sunlight. Do not expose to temperatures in excess of 50°C/122°F.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.

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